



# ELECTROLAST™

At Home Micro-Current Facial Device

Minimizes facial lines and folds  
Reduces under eye puffiness  
Lifts and firms the skin  
Improves skin tone

Operating Manual

# INTRODUCTION

Cosmetic science reveals that micro-current facials have positive anti-aging benefits. This method is a holistic alternative to improve the appearance of the skin and is safer and gentler than invasive procedures like plastic surgery or controversial injectable substances such as Botox. I have performed micro current facials on so many happy customers in my spa practice and found there is an improvement to overall-muscle and skin tone. In response to my own customers asking for a home version, I have created ElectroLast. Now everyone can have the "Red-Carpet" micro-current treatment to maintain the health of the skin. Best of all, ElectroLast is easy to use, travels with you and is affordable.

The ElectroLast home device adopts low electronic, gentle current that helps stimulate circulation, minimize fine lines and lift sagging skin. A slight tingling sensation is felt and lets you know the machine is working on muscles. The intensity output is adjustable so at no point should you feel uncomfortable. You are in control to adjust the frequency to your personal comfort level. Overtime as your skin adjusts to the device the sensation may diminish and sometimes you may not feel any sensation at all. Please note however that the machine is still working on stimulating the facial muscles.

Soon you will see the effects of ElectroLast. Your skin will look firmer, your facial muscles toner, reduced lines and skin folds, and an overall glow to your complexion.

Join me as an "age-rebel" and enjoy your youthful appearance.

*Best, Marlene*

*Marlene Katz is an esthetician and spokesperson for the cosmetic industry. Her beauty tips and articles have been featured on Good Morning America, Better TV, Elle Magazine, Harpers Bazaar and many trade publications.*



## Read all instructions before using this appliance

Use this product only as directed in this manual. This product is not intended for children. This product is not to be used inside the eyes, on the eye lid or on the eyelashes. It is made for the face and not to be used on the body. Do not use when driving. If you are under medical care, please check with your physician for approval before using this product.

### Do not use:

- If you have an acute disease or are epileptic
- When pregnant, undergoing chemotherapy or radiation treatments
- If you have a pace-maker, metal prosthesis or facial implants.
- If you have skin trauma such as open wounds, dermatitis or skin cancer.
- If you have inflammation or open skin wounds herpes simplex on the face or neck, skin rash, acne eruptions, sunburn
- If you are on Accutane
- If you are recovering from Plastic Facial Surgery
- If you have skin bruising
- If you have had any of the following in the past few weeks: laser therapy, chemical peels, reddening of the skin from Retin-A applications or peels.

**Mode I: Select and Red LED light flashes**

Current Consumption 3.5mA to 4.7mA

Output Current: 1uA-2uA

Voltage: 18V to 22V

Output Frequency: 30ms-3333ms (0.3 to 33Hz)

**Mode II: Select and Blue LED light flashes**

Current Consumption: 3.5 mA to 4.7 mA

Output Current: 1uA-50uA

Voltage: 20V to 22V

Output Frequency: 11ms (90 Hz)

**Mode III: Select and Violet light flashes**

(Blue and Red LED light flashing together)

Current Consumption: 3.5mA to 4.8 mA)

Output Current: 1uA-30uA

Voltage: 20V to 35V

Output Frequency: 36ms (28 Hz)


Output Frequency for LED: five levels of intensity

The three modes have the same level of output of frequency for the LED levels of intensity.

1. 0.488 Hz / 2. 0.606 Hz / 3. 0.8 Hz / 4. 1.176 Hz / 5. 2.22 Hz

(\*Note: uA is 'micro-amp and 1ms=1000us=1000\*1000ns=1000/Hertz)


Size: 30(dia) 18 (H) mm Weight: about 14g

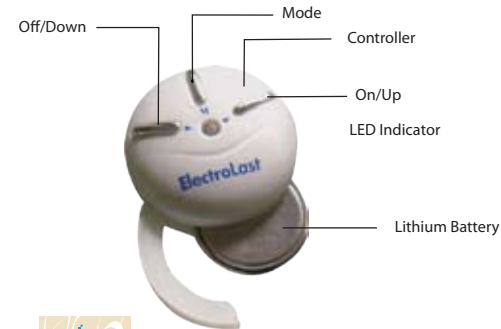

Manual



Pouch

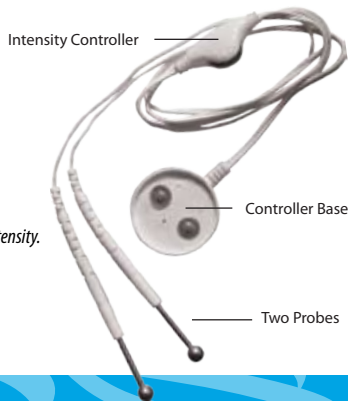

CE RoHS Certified

## Product Specifications

To Begin Operations with ElectroLast

1. Affix the Controller to the circular Controller Base. The controller easily snaps into the magnetic base.
2. Have the controller facing you with the three operating buttons on top.
3. Push the button to the right one time to turn on the device. When the power is on, the red LED indicator will light. The mode is set to level 1. There are five levels of intensity with level 1 as the lowest of frequency and level 5 as the highest.

*Note: If you turn on the device but do not begin use, it will automatically shut off in one minute.*

4. To increase the intensity press down on the right button. The level begins at 1 and increases to 5 with each press. To decrease the intensity, press the left button down. Each press on the left button reduces the intensity down. (If you are on level five then press once to reduce to level four and so on. The sixth press will turn the unit off).
5. To turn the unit off at any level, press the left button down until the level is at 1 and then once more to turn the device off. When the device is turned off the LED indicator light turns off.
6. The middle button controls the LED light modes. There are three LED light modes.

### Mode I - Red LED Light

Red indicator light mode is used to strengthen the longer facial muscles.

### Mode II - Blue LED Light

Blue indicator is used to strengthen the shorter facial muscles.

### Mode III - Violet LED Light

Violet indicator is used to stimulate the circulation and minimize fine facial lines and folds.

## FACIAL INSTRUCTIONS

1. Cleanse the face and eye area thoroughly with a mild face cleanser and rinse well with water.
2. The micro-current works best on damp skin. To keep the skin wet, have a bowl of water and a cotton pad or sponge set aside to wet the face as needed.
3. Apply Lock-In Cap™ to face and neck. Use one to three drops. Massage into face and neck.

*Suggested Levels of Intensity: Please note that the suggested levels can be adjusted as needed to suit your own comfort level. At no time, should the treatment be uncomfortable.*

### Mode I - Red LED

Sliding: These movements address the longer muscles. Remember to begin at the base of the neck and slide the probes upwards.

~ Use level 2 to 3 for first pass around face.

### Mode II - Blue LED

Lift & Hold: These movements address the shorter muscles.

Gently lift the skin up with each pinching movement. Note: The diagram uses lines to show whether to hold the machine vertical or to hold it horizontal.

~ Use level 3-4 for second pass around face.

### Mode III - Red & Blue LED

Erase & Feather: Use Erase & Feather strokes to smooth away fine lines.

Zig Zag  Criss Cross  XXXX  Circles 

Use short and light feather type movements like a brush stroke with water colors.

~ Use level 4-5 for the third pass around face.

## EYE INSTRUCTIONS

~ Use level 1-2 for all passes around the eye.

*\*Caution: Do not put product and do not use the ElectroLast™ device on or close to tear ducts or on top of upper eye lids.*

## Facial Movements

Always begin at the base of the neck. Work on one side of the face at a time and follow the diagram. You are always working against gravity moving the probes upwards and/or outwards, never down.

Please follow the diagram for proper usage.

Step 1. Base of neck to jaw line

Step 2. Under chin to jaw

Step 3. Chin to lower lip line

Step 4. Lower jaw to lower cheek


Step 5. Upper lip area

Step 6. Lower cheek bone to upper cheek bone

Step 7. Lower outer nostril to center of brow

Step 8. Outside nose area to eye areas

Step 9. Forehead to hair

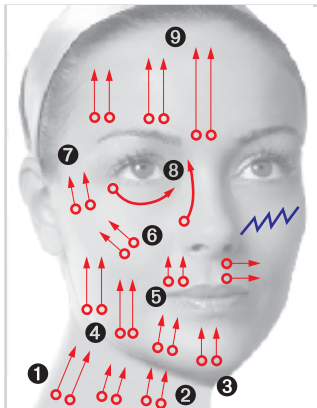


SLIDE



2 wand tips



1 wand tip

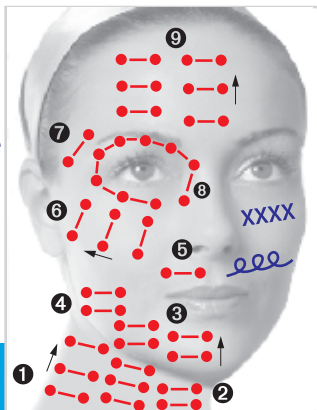

LIFT & HOLD



2 wand tips



1 wand tip



direction

## Frequently Asked Questions

*Q. How often can I use the machine?*

*A. You can use it as often as three times a week with a day off in between.*

*Example: Use ElectroLast Monday, Wednesday, and Saturday.*

*Q. Why can't I use it everyday?*

*A. When the muscles rest they work best. Over exercising the muscle is not recommended.*

*Q. What if I do not feel the current?*

*A. You may need to re-mist the skin as the current works best when the skin is wet. You can increase the current to one higher level on the face. Eyes must always stay on #1, 2 or 3 for all passes.*

*Q. What if I feel a muscle twitch?*

*A. It is normal to feel a slight tingle from the current and a slight warming of the skin. However, if the muscle twitches, reduce the current. There are reasons that the skin may be sensitive such as hormones, weather changes and stress. Adjust the current up or down so it is comfortable for you.*

*Q. Should my skin be red after the treatment?*

*A. A rosy glow is expected as the ElectroLast stimulates the circulation. Reddened, rashes or irritated skin are not a normal reaction and if you experience this, stop using the machine and contact Customer Service: [info@intelligentskinsense.com](mailto:info@intelligentskinsense.com) or 800-544-0040.*

### **Maintenance of Unit:**

1. Wipe the two probes tips clean after each use with dish soap and water. Note: Do not ever immerse the controller or the wires in water.
2. Dust the controller clean with a tissue after each use. If need be, use alcohol on the tissue to clean the device.
3. Keep product contained in the travel pouch or in the box when not in use. Do not store in high temperatures.

Warranty: The product is warranted for one year from the purchase date and only with a receipt to verify date of purchase. If the machine needs repair, contact customer service at 800-544-0040 or [info@intelligentskinsense.com](mailto:info@intelligentskinsense.com).